

	Jet Compact / Jet Wire	Home Jet	AMBUJET
Dimensions	245 x 70 x 185 mm	200 x 225x 285 mm	200 x 250x 320 mm (with fixing)
Weight	1.6 kg	2.4 kg	2.4 kg + 0.8 kg (fixing)
Vacuum range	from 50 to 550 mmHg	from 50 to 550 mmHg	from 50 to 550 mmHg
Flow (free flow)	32 l/min	32 l/min	32 l/min
Noise level	< 70db	< 70db	< 70db
Electrical connection requirements	12 VDC 33 W max	12 VDC 33W max	12 VDC 33W max
Battery	12 V DC Ni-MH	12 V DC Ni-MH	12 V DC Ni-MH
Battery operating time	max 45-60 min at free flow	max 45-60 min at free flow	max 45-60 min at free flow
Operating time	30 min. ON, 30 min. OFF	30 min. ON, 30 min. OFF	30 min. ON, 30 min. OFF
Battery recharging time	5 hours	5 hours	5 hours
ISO Classification	High vacuum/High flow	High vacuum/High flow	High vacuum/High flow
IP protection	IP12	IP12	IP12
Medical device classification	Ila	Ila	Ila
Patient tube	1.0 m	1.8 m	1.8 m
Operating temperature	From 0°C to +45°C	From 0°C to +45°C	From 0°C to +45°C
Storage temperature	From - 40°C to +70°C	From - 40°C to +70°C	From - 40°C to +70°C
Umidity	0 to 95%	0 to 95%	0 to 95%
Charging options	Power plug for 12 V	Power plug for 12 V	Rechargeable fixing
Other options	NATO stock number coded	Canister protection from 1.2 m falls	20 G UNI EN-1789 Crash Test certificate
Canister capacity	300 ml disposable, 1000 ml reusable, 1000 ml Serres	800 ml disposable, 1000 ml reusable, 1000 ml Serres	800 ml disposable, 1000 ml reusable, 1000 ml Serres

UNBEATABLE FEATURES ON THE MARKET

Find out which **Spencer suction device** best suits your needs

Model	Canister type	Shock-proof	Rechargeable fixing	10 G fixing	20 G fixing
HomeJet SC75200D	800 ml	•			
HomeJet SC75300D	1000 ml	•			
HomeJet SC75210B	Serres 1000 ml	•			
Jet Compact SC75100D	300 ml				
Jet Wire SC75410C	Serres 1000 ml				
Jet Wire SC75400C	1000				
AMBUJET SC75560D	Serres 1000 ml	•	•	•	
AMBUJET SC75510D	Serres 1000 ml	•	•		•
AMBUJET SC75760B	Serres 1000 ml	•		•	
AMBUJET SC75710B	Serres 1000 ml	•			•
AMBUJET SC75550C	800 ml	•	•	•	
AMBUJET SC75500C	800 ml	•	•		•
AMBUJET SC75750C	800 ml	•		•	
AMBUJET SC75700C	800 ml	•			•
AMBUJET SC75650D	1000 ml	•	•	•	
AMBUJET SC75600C	1000 ml	•	•		•
AMBUJET SC75850C	1000 ml	•	•	•	
AMBUJET SC75800C	1000 ml	•	•		•
AMBUJET S SC75870C	1000 ml	•		•	
AMBUJET S SC75770C	Serres 1000 ml	•		•	
AMBUJET S SC75670C	1000 ml	•	•	•	
AMBUJET S SC75570C	Serres 1000 ml	•	•	•	


AMBUJET 20 G FASTENERS

Easy to install and one-click operation

Downtime is the enemy of productivity. That is why AMBUJET is engineered to be easy and fast to install, while minimizing the need for routine maintenance. Thanks to the perfect integration between the fastener and the suction unit, the overall dimensions remain almost identical to those of the suction unit itself.



Quick release button

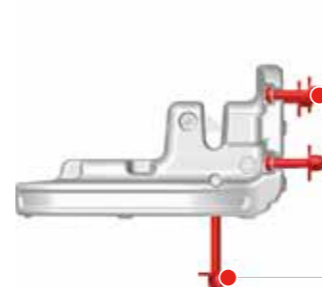

Two types of fixing


Wall mount

Mount for surface

AMBUJET is the only suction unit on the market which can be connected to the base and operated with one-hand, making such device the best choice for ambulance and special vehicles builders. **AMBUJET is the first electrical suction machine to obtain 20 G UNI EN-1789 Crash Test certificate.** A record that underline our attention towards details and the high technological value of the design.

DISCOVER THE NEW AMBUJET S WITH 10 G FASTENER

Reliable and interchangeable

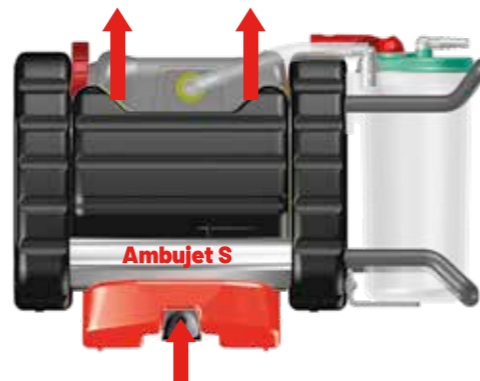

Quick release button



Quick release button

Innovation starts with understanding what our customers need and, as a global leader in EMS technologies, **Spencer is anticipating all the challenges of a changing industry.** Thanks to our latest innovation, we create a **new layout for fasteners** which are fully interchangeable between **Ambujet S** and **all NXT ventilators** made by Spencer.

The value of a choice is measured in the future.

Jet series

Introducing the ultimate expression of Spencer's efficiency, design and innovation.

spencer.it



SPENCER®

JET COMPACT

Portable compact suction device

- COMPACT
- PRACTICAL
- LIGHTWEIGHT



Jet Compact is ideal for use in the emergency field, during transport or in the hospitals. You will be surprised by its compact and lightweight versatility. The device is also appreciated in pre-hospital and homecare sector for its low level of noise and vibrations.

Unique features

- Only 1.6 kg with high performance
- Compact and easy to fit in an emergency bag
- DO-160 airborne equipment compatibility and NATO stock number coded

Standard equipment

- 220 V battery
- Single canister and antibacterial filter
- Disposable patient tube
- 12 V cigarette lighter cable
- Shoulder strap

Optional accessories

- AC/DC adapter 100/240 VAC 50-60 Hz
- Transport Bag

LED display for use in night or poor weather conditions

High performance with intuitive and easy to use controls

High portability which allows to be positioned in a rucksack or in an emergency bag

No tools required to change the battery



HOMEJET

Portable suction unit

- SHOCK PROOF
- INDESTRUCTIBLE
- STABLE



HomeJet is able to offer a suction capacity that is unmatched in its class. It's the only suction device with fast battery substitution with no tool required. It allows for easy and quick replace of the canister that is directly connected to the suction unit so as to reduce the number of components used and also any repair costs.

Unique features

- Only 2.4 kg with full canister protection
- Rapid substitution of the battery
- Direct connection between canister and suction device

Standard equipment

- 220 V battery
- Single canister and antibacterial filter
- Disposable patient tube
- 12 V cigarette lighter cable
- Shoulder strap

Optional accessories

- AC/DC adapter 100/240 VAC 50-60 Hz
- Transport Bag

High suction capacity of about 32 liters/minute

Easy conversion between different canisters

Low level of noise and vibration



Polyethylene shell to withstand any impact

AMBUJET

Portable suction unit

- 20 G TESTED
- QUICK ASSEMBLY
- CANISTER PROTECTION



Durable, lightweight, stable, compact. In just one word AmbuJet. The only device with an integral indestructible polyethylene shell able to fully protect both the unit and the canister from impact damage, while its wide base makes the device exceptionally stable. AmbuJet is the first suction unit equipped with a certified 20 G fixing.

Unique features

- Wide base for increased stability
- Time-saving thanks to easy and quick installation
- 20 G crash test certificate

Standard equipment

- Fixing
- 220 V Battery
- Single canister and antibacterial filter
- Disposable patient tube
- Fixing with integrated charging system
- Shoulder strap

Optional accessories

- AC/DC adapter 100/240 VAC 50-60 Hz
- Transport Bag

Direct device canister connection to facilitates insertion and substitution of the canister

Easy to install and one-click operation

Wide base for **exceptional stability**



Small and lightweight with **one-hand-grip design**

IDEAL CHOICE FOR EMS

IDEAL CHOICE FOR HOMECARE

IDEAL CHOICE FOR AMBULANCE